

# Oxford<sup>®</sup>



## Presence



 **Joerns**  
Healthcare

Redefining patient handling

Outstanding performance  
and lifting range

# Presence

'Presence' is a lifting device that is competent enough to handle almost any patient-handling task. It is an all-round performer that gives you more than you expect!

## Key Features:

- Outstanding lifting range (1910mm)
- Greater spatial area for the carer and patient
- Powered base
- Large lifting capacity (SWL) 227kgs/35st
- Interchangeable spreader bar/cradle configuration
- Powered cradle (optional)
- Digital weigh-scale (optional)

In addition to this versatility, the Presence has been developed using Oxford's key design principles of;

- Ergonomic considerations
- Intelligent positioning of key interaction points

These combined elements mean the Presence lift meets the true needs of the environment, carer and patient.



Convenient hand control rest



Ergonomic foot push pad



Ergonomic oversized push handle

## Outstanding lifting range

Whether from the floor or onto a high surface the Presence has an astonishing range of movement. The minimum height with the six point spreader bar is just 440mm, which ensures excellent off the floor transfers.



Conversely, if you wish to access high beds or nursing surfaces the Presence will lift up to 1910mm. This is especially useful where you have variable height beds, special lift mattresses and other high surfaces. This type of flexibility is increasingly important in the care environment. It is the difference between effective safe transfers and unnecessary risk for both the carer and patient.



## Greater spatial area for carer and patient



Sculptured base design allows carer to get close to the lift, which promotes ease of use. The additional space promotes better ergonomics for the carer, as it is possible to interact more closely with the patient.



## Powered base

The Presence has two powerful electric motors which open and close the legs automatically. Additionally, the Presence's unique 'swan-neck' leg design also allows the lift to get close to the widest of obstacles. This is especially useful when encountering large chairs, bulky commodes and wheelchairs.



### **Tapered 'swan-neck' leg design**

Allows greater internal area, whilst not compromising overall product width.



## Large lift capacity

In today's demanding care environment, larger patients are becoming more of a challenge for conventional lifts.


The Presence provides a larger spatial area to accommodate these patients. Transferring now becomes safer and easier to achieve.

Plus, with a mighty safe working load of 227kg/35st, the Presence lift is able to cope with the vast majority of patient handling situations.

## Interchangeable spreader bar/4 point positioning cradle with weigh-scale options



(A)  
6-point spreader bar

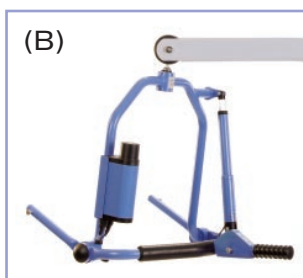


(B)  
Manual or powered cradles



Digital weigh-scale  
(Class III)

One of the restrictions of using certain lifts is the type of slings you can use on the product. The Presence, however, gives you a choice of sling systems because you can use the conventional 6-point spreader bar (A) or Oxford 4-point positioning cradle (B - powered version shown).

This provides the carer with a greater choice of slings when considering support, comfort and flexibility for the patient. The 4-point cradle is also available in both manual and powered versions.

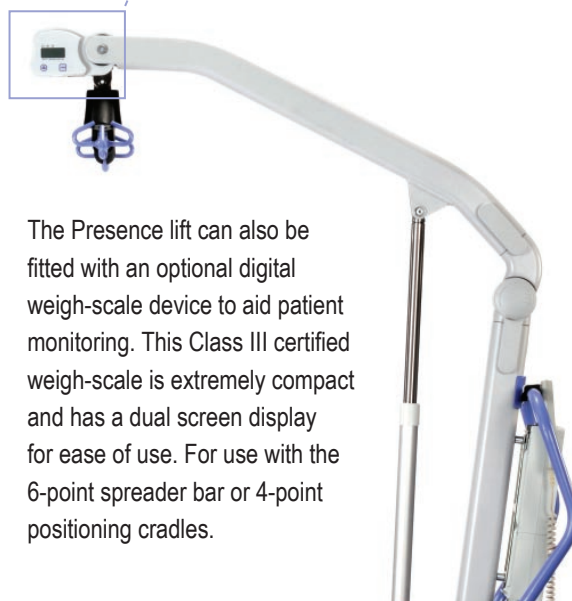

The Presence lift can also be fitted with an optional digital weigh-scale device to aid patient monitoring. This Class III certified weigh-scale is extremely compact and has a dual screen display for ease of use. For use with the 6-point spreader bar or 4-point positioning cradles.

Specification	Imperial	Metric
Safe Working Load	500 lbs	227 kgs
Maximum Overall Length	57.5 inches	1460 mm
Minimum Overall Length	55.1 inches	1400 mm
Maximum Overall Height	81.3 inches	2065 mm
Minimum Overall Height	57.5 inches	1460 mm
Spreader Bar - Maximum Height	75.2 inches	1910 mm
Spreader Bar - Minimum Height	16.9 inches	430 mm
Height at Maximum Reach	51.5 inches	1310 mm
Reach at Maximum Height	25.6 inches	650 mm
Reach at Minimum Height	15.7 inches	400 mm
Maximum Reach (centre of spreader bar to front of the mast)	36.6 inches	930 mm
Turning Radius	64.2 inches	1630 mm
Legs Open - External Width	45.2 inches	1150 mm
Legs Open - Internal Width	40.1 inches	1020 mm
Legs Closed - External Width	26.4 inches	670 mm
Legs Closed - Internal Width	21.6 inches	550 mm
Overall Height of Legs	4.7 inches	120 mm
Ground Clearance	1.4 inches	35 mm
Front Twin Castors	4.0 inches	100 mm
Rear Braked Castors	4.0 inches	100 mm
4 Point Positioning Cradle - Maximum Height	61.0 inches	1550 mm
4 Point Positioning Cradle - Minimum Height (usable)	15.4 inches	390 mm

Weights	Imperial	Metric
Mast, Base & Boom Assembly (includes 6 point spreader bar)	81.4 lbs	37.0 kgs
Power Pack	6.6 lbs	3.0 kgs
Total (including 6-point spreader bar)	88.0 lbs	40.0 kgs
Base Assembly (excluding battery)	44.0 lbs	20.0 kgs
Mast & Boom Assembly (excluding battery)	37.4 lbs	17.0 kg
4 Point Positioning Cradle (manual operation)	13.2 lbs	6.0 kgs
4 Point Positioning Cradle (powered operation)	19.8 lbs	9.0 kgs

Specifications may vary due to manufacturing tolerances

#### Electrical Specifications (Lift)

Battery 1 x 12 volt rechargeable sealed lead acid type  
 Battery capacity 3.2 Ah  
 Charger rated input 230V AC 50/60Hz  
 Charger rated output 27.4/29.0 V DC@0.8A

#### Electrical Specifications (Powered cradle only)

Battery 24V DC rechargeable NiMH sealed  
 Battery capacity 1400 mAh  
 Charger rated input 100 - 240V AC / 24 V DC 50/60Hz Max. 75mA  
 Charger rated output Max 28.4V DC

#### Warranty

The Presence comes with a 2 year warranty as standard.  
 Please refer to your user manual for further information.

#### Electrical Shock Protection

Charger - class II  
 Lift - internal power source

Degree of shock protection  
 Charger - Type B  
 Lift - Type B  
 Intended operating environment : >+5° <+40°  
 Outside this environment functionality and safety may be compromised



# Oxford®



Joerns Healthcare Limited  
 Drakes Broughton Business Park, Worcester Road, Drakes Broughton  
 Pershore, Worcestershire, WR10 2AG United Kingdom  
 (T) +44 (0) 844 811 1156 (F) +44 (0) 844 811 1157  
 www.joerns.co.uk • email: info@joerns.co.uk