

# Oxford®



## Stature



 **Joerns**  
Healthcare

Redefining patient handling

Outstanding all-round performer

# Stature

The 'Stature' patient lift, Oxford's flagship product, recognises the true needs of the modern care environment.

## Key Features:

The Stature has many key design features which make it a truly flexible piece of equipment:

- Large lift capacity (SWL) 227kgs/35st
- Outstanding lifting range (Min 390mm, max 1550mm)
- Interchangeable spreader bar/cradle configuration
- Powered base
- Powered cradle (optional)
- Digital weigh-scale (optional)

## Outstanding lifting range



The Stature's incredible lifting range enables it to cope easily with the most critical of daily tasks:

- Chair to chair transfers
- Bed to chair transfers
- Floor to chair transfers
- Chair to bath transfers

## Interchangeable spreader bar/positioning cradle




Like all Oxford Pro passive mobile lifts, you have the option to choose whichever cradle or spreader bar best suits the needs of the care environment and patient alike. With a wide selection of slings and attachments available, you can be sure to get the right combination each and every time.





The Stature comes complete with a 4-point positioning cradle as standard. This enables the carer to position the patient perfectly, simply by rotating the cradle (and patient) into a fully seated or reclined orientation.

The 4-point positioning cradle is available in manual or powered versions. The powered cradle allows the carer to lower and position the patient independently for optimum recline and comfort - eliminating the requirement for manual repositioning.

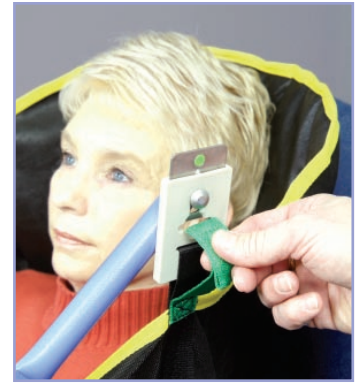

**Sling safety:** With Oxford's unique 'Securi3' sling attachment system you can be assured that patient safety is not compromised. The Securi3 system ensures that no inadvertent detachment occurs, once the sling clip has been attached to the cradle. This guarantees that each transfer will be achieved safely and efficiently each and every time.

## Powered base

Opening and closing of the legs is fully automated and leaves the carer free to focus on the requirements of the patient.



## Extensive range of accessories

The Stature is available with a number of key accessories, making it both versatile and flexible in use.



6-point spreader bar



Manual or powered cradles



Digital weigh-scale (Class III certified)

**5** The Stature was designed with the 5 key principles of moving and handling in mind. This direct correlation ensures Oxford's lift designs are simple, safe and above all comfortable to use.

This view is shared by our clinical consultant teams, who have assisted us in conducting detailed usability surveys to ensure that the Stature meets the true needs of the patient, the carer and the environment in which it operates.

- 1 Get close to the load
- 2 Use a wide stable base
- 3 Ensure a comfortable firm grasp
- 4 Keep spine close to neutral
- 5 Make sure movement is smooth

Specification	Imperial	Metric
Safe Working Load	500 lbs	227 kgs
Maximum Overall Length	54.7 inches	1390 mm
Minimum Overall Length	50.8 inches	1290 mm
Maximum Overall Height	83.8 inches	2130 mm
Minimum Overall Height	76.6 inches	1945 mm
Maximum Height to Attachment Point - 4 Point Cradle	61.0 inches	1550 mm
Minimum Height to Attachment Point - 4 Point Cradle	15.4 inches	390 mm
Maximum Height to Attachment Point - 6 Point Spreader Bar	77.7 inches	1975 mm
Minimum Height to Attachment Point - 6 Point Spreader Bar	26.6 inches	675 mm
Turning Radius	61.0 inches	1550 mm
Legs Open - External Width	53.5 inches	1360 mm
Legs Open - Internal Width	48.4 inches	1230 mm
Legs Closed - External Width	29.9 inches	760 mm
Legs Closed - Internal Width	25.1 inches	640 mm
Overall Height of Legs	4.3 inches	110 mm
Ground Clearance	1.2 inches	30 mm
Front Twin Castors	4.0 inches	100 mm
Rear Braked Castors	4.0 inches	100 mm

Weights	Imperial	Metric
Mast, Base & Boom Assembly	142 lbs	64.5 kgs
4 Point Positioning Cradle (manual operation)	13.2 lbs	6.0 kgs
Power Pack	6.6 lbs	3.0 kgs
Total	162 lbs	73.5 kgs
4 Point Positioning Cradle (powered operation)	19.8 lbs	9.0 kgs

Specifications may vary due to manufacturing tolerances

#### Electrical Specifications (Lift)

Battery 1 X 12 volt rechargeable sealed lead acid type  
 Battery capacity 3.2 Ampere hours  
 Charger rated input 230Vac 50/60Hz  
 Charger rated output 27.4/29.0 VDC@0.8A

#### Electrical Specifications (Powered cradle only)

Battery 24V DC rechargeable NiMH sealed  
 Battery capacity 1400 mAh  
 Charger rated input 100 - 240V AC / 24 V DC 50/60Hz Max. 75mA  
 Charger rated output Max 28.4V DC

#### Warranty

The Stature comes with a 2 year warranty as standard.  
 Please refer to your user manual for further information.

#### Electrical Shock Protection

Charger - class II  
 Lift - internal power source

Degree of shock protection  
 Charger - Type B  
 Lift - Type B

Intended operating environment : >+5° <+40°  
 Outside this environment functionality and safety may be compromised



# Oxford®



Joerns Healthcare Limited  
 Drakes Broughton Business Park, Worcester Road, Drakes Broughton  
 Pershore, Worcestershire, WR10 2AG United Kingdom  
 (T) +44 (0) 844 811 1156 (F) +44 (0) 844 811 1157  
 www.joerns.co.uk • email: info@joerns.co.uk